

Initial Wi-Fi & Credentials Setup

FIRST-TIME USE

1 Power On

Plug in your SugrBee monitor.

2 Connect to Wi-Fi

On your phone or PC, connect to "SugrBee Setup" (no password required). The portal should open automatically. If not, toggle Wi-Fi off/on or enter the IP shown on screen (e.g., 192.168.4.1).

3 Configure Credentials & Settings

Tap Configuration. Select home Wi-Fi and enter its password. Enter Dexcom or Libre credentials, then customize colors and alert thresholds.

4 Save & Restart

Tap Save. Device restarts and shows glucose data.

FOR LIBRE USERS ONLY

Libre users must be following the account providing the credentials in SugrBee. You must follow yourself if you are using your account.


1. Open the primary FreeStyle Libre app.
2. Tap Profile (bottom right) or the top menu.
3. Select Connected Apps → LibreLinkUp.
4. Tap Add Connection and invite your own email address.
5. Accept the invite and complete signup.
6. Enter LibreLinkUp credentials into SugrBee setup.



Touchscreen Controls

Dot Icon
Opens the main Settings Menu.

Glucose Number
Switches unit between mg/dL and mmol/L.

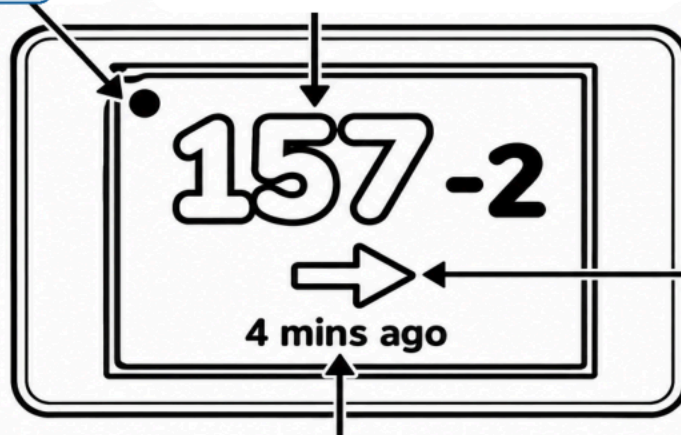

Trend Arrow
Toggles between trend arrow and 2 inline graph views.

Time Display
Switches between time since last reading and current time/date.

Settings Menu Overview

Tap the top-left dot to open Settings; tap the top-right dot inside Settings for Page 2.

SETTING	WHAT IT DOES
Brightness	Increases or decreases screen brightness intensity.
Night Mode	Dims screen to lowest possible brightness for sleeping.
Invert Colors	Swaps background and text colors for dark/light rooms.
Rotate Display	Flips display 180° for convenient power cable positioning.
Update	Checks for updates; button turns green when available.
Setup	Re-opens initial setup without erasing saved credentials.
Reset	Clears saved settings and returns to first-time setup.
Down Detector	Checks Dexcom or Libre server status for outages.
Debug	Sends a diagnostic log to customer support.

Troubleshooting Tips

ISSUE	WHAT TO TRY
No Glucose Data	· Verify sensor is active and transmitting. · "No Data" appears if offline for 6+ minutes. · Check Wi-Fi/server status via Down Detector. · Unplug for 10 seconds and reconnect (or use Reset).
Incorrect Time	· Device automatically syncs time based on location. · Unplug for 30 seconds and reconnect to re-sync.
Can't Connect to "SugrBee Setup"	· Network may take 2-5 minutes to appear; toggle Wi-Fi. · If portal does not open, enter the displayed IP. · Disable "Smart Wi-Fi Switch" on Android.
Can't Log In	· Dexcom: Enable "Dexcom Share" in the Dexcom app. · Libre: Set up LibreLinkUp and follow yourself.

COMMUNITY & SUPPORT

Facebook: SugrBee Glucose Monitor (Official) | YouTube: Type1Treasures | Instagram: @Type1Treasures